

DTIC COPY



AD-A229 952

TECHNICAL REPORT  
NATICK/TR-91/004

AD \_\_\_\_\_

# A COMPUTER PROGRAM FOR MEASURING FIBERS WITH THE ZEISS CSM 950 SCANNING ELECTRON MICROSCOPE

Deborah Sakelakos  
Samuel H. Cohen

November 1990

Final Report  
June 1990 to September 1990

DTIC  
ELECTE  
DEC 11 1990  
S E D

Approved for public release, distribution unlimited

US ARMY NATICK RESEARCH, DEVELOPMENT AND  
ENGINEERING CENTER  
NATICK, MA 01760-5000

SOLDIER SCIENCE DIRECTORATE

00 10 11 1990

## DISCLAIMERS

The findings contained in this report are not to be construed as an official Department of the Army position unless so designated by other authorized documents.

Citation of trade names in this report does not constitute an official endorsement or approval of the use of such items.

## DESTRUCTION NOTICE

For Classified Documents:

Follow the procedures in DoD 5200.22-M, Industrial Security Manual, Section II-19 or DoD 5200.1-R, Information Security Program Regulation, Chapter IX.

For Unclassified/Limited Distribution Documents:

Destroy by any method that prevents disclosure of contents or reconstruction of the document.

<b>REPORT DOCUMENTATION PAGE</b>		Form Approved OMB No 0704-0188	
<small>1. This report contains information which has been determined by the originating agency to be exempt from automatic declassification under Executive Order 11652. It is therefore exempt from automatic declassification under E.O. 11652. The exemption is authorized by 5 U.S.C. 552(a)(1) and (2). The exemption is subject to review by the Department of Defense Information Security Management System (DODISMS) and the Department of Defense Information Security Management System (DODISMS) Report of Exemption (DOE) Form 10704-0188 Washington, DC 20503.</small>			
<b>1. AGENCY USE ONLY (leave blank)</b>	<b>2. REPORT DATE</b> 5 November 1990	<b>3. REPORT TYPE AND DATES COVERED</b> Final Report - June 1990 to Sep 1990	
<b>4. TITLE AND SUBTITLE</b> A Computer Program for Measuring Fibers with the Zeiss CSM 950 Scanning Electron Microscope		<b>5. FUNDING NUMBERS</b>  PE: 612786 PR: 1L162786AH99BAC00 WU: AH99BAC00	
<b>6. AUTHOR(S)</b> Deborah Sakelakos and Samuel H. Cohen			
<b>7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES)</b> U.S. Army Natick Research, Development and Engineering Ctr Kansas Street ATTN: STRNC-YSP Natick, Massachusetts 01760-5020		<b>8. PERFORMING ORGANIZATION REPORT NUMBER</b>  NATICK/TR-91/004	
<b>9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES)</b>		<b>10. SPONSORING/MONITORING AGENCY REPORT NUMBER</b>	
<b>11. SUPPLEMENTARY NOTES</b> Deborah Sakelakos was a participant in the FY90 U.S. Army High School and Mathematics Faculty Program funded by the Army Research Office.			
<b>12a. DISTRIBUTION AVAILABILITY STATEMENT</b>  Approved for Public Release, Distribution Unlimited		<b>12b. DISTRIBUTION CODE</b>	
<b>13. ABSTRACT (Include keywords)</b>  A computer program designed for the stereological measurement of fibers was written for use with a ZEISS CSM 950 Scanning Electron Microscope (SEM). This program utilized the ZEISS' build-in Kontron Image Analyzer to access measurements of a variety of parameters for fiber analysis of images from a video camera, light microscope or SEM.  P. J.			
<b>14. SUBJECT TERMS</b> Computer program, ZEISS CSM 950, Kontron Image Analyzer Scanning Electron Microscope, Fibers, Images, Parameters		<b>15. NUMBER OF PAGES</b> 26 <b>16. PRICE CODE</b>	
<b>17. SECURITY CLASSIFICATION OF REPORT</b> UNCLASSIFIED	<b>18. SECURITY CLASSIFICATION OF THIS PAGE</b> UNCLASSIFIED	<b>19. SECURITY CLASSIFICATION OF ABSTRACT</b> UNCLASSIFIED	<b>20. LIMITATION OF ABSTRACT</b> UNLIMITED

## GENERAL INSTRUCTIONS FOR COMPLETING SF 298

The Report Documentation Page (RDP) is used in announcing and cataloging reports. It is important that this information be consistent with the rest of the report, particularly the cover and title page. Instructions for filling in each block of the form follow. It is important to *stay within the lines* to meet optical scanning requirements.

**Block 1. Agency Use Only (Leave blank).**

**Block 2. Report Date.** Full publication date including day, month, and year, if available (e.g. 1 Jan 88). Must cite at least the year.

**Block 3. Type of Report and Dates Covered.** State whether report is interim, final, etc. If applicable, enter inclusive report dates (e.g. 10 Jun 87 - 30 Jun 88).

**Block 4. Title and Subtitle.** A title is taken from the part of the report that provides the most meaningful and complete information. When a report is prepared in more than one volume, repeat the primary title, add volume number, and include subtitle for the specific volume. On classified documents enter the title classification in parentheses.

**Block 5. Funding Numbers.** To include contract and grant numbers; may include program element number(s), project number(s), task number(s), and work unit number(s). Use the following labels:

C - Contract	PR - Project
G - Grant	TA - Task
PE - Program Element	WU - Work Unit Accession No.

**Block 6. Author(s).** Name(s) of person(s) responsible for writing the report, performing the research, or credited with the content of the report. If editor or compiler, this should follow the name(s).

**Block 7. Performing Organization Name(s) and Address(es).** Self-explanatory.

**Block 8. Performing Organization Report Number.** Enter the unique alphanumeric report number(s) assigned by the organization performing the report.

**Block 9. Sponsoring/Monitoring Agency Name(s) and Address(es).** Self-explanatory.

**Block 10. Sponsoring/Monitoring Agency Report Number.** (If known)

**Block 11. Supplementary Notes.** Enter information not included elsewhere such as: Prepared in cooperation with...; Trans. of...; To be published in.... When a report is revised, include a statement whether the new report supersedes or supplements the older report.

**Block 12a. Distribution/Availability Statement.** Denotes public availability or limitations. Cite any availability to the public. Enter additional limitations or special markings in all capitals (e.g. NOFORN, REL, ITAR).

DOD - See DoDD 5230.24, "Distribution Statements on Technical Documents."

DOE - See authorities.

NASA - See Handbook NHB 2200.2.

NTIS - Leave blank.

**Block 12b. Distribution Code.**

DOD - Leave blank.

DOE - Enter DOE distribution categories from the Standard Distribution for Unclassified Scientific and Technical Reports.

NASA - Leave blank.

NTIS - Leave blank.

**Block 13. Abstract.** Include a brief (*Maximum 200 words*) factual summary of the most significant information contained in the report.

**Block 14. Subject Terms.** Keywords or phrases identifying major subjects in the report.

**Block 15. Number of Pages.** Enter the total number of pages.

**Block 16. Price Code.** Enter appropriate price code (*NTIS only*).

**Blocks 17. - 19. Security Classifications.** Self-explanatory. Enter U.S. Security Classification in accordance with U.S. Security Regulations (i.e., UNCLASSIFIED). If form contains classified information, stamp classification on the top and bottom of the page.

**Block 20. Limitation of Abstract.** This block must be completed to assign a limitation to the abstract. Enter either UL (unlimited) or SAR (same as report). An entry in this block is necessary if the abstract is to be limited. If blank, the abstract is assumed to be unlimited.

# TABLE OF CONTENTS

	<u>Page</u>
LIST OF FIGURES	iv
LIST OF TABLES	v
PREFACE	vii
INTRODUCTION	1
METHODS	1
RESULTS	12
CONCLUSIONS	12
REFERENCE	20

Accession For	
NTIS GRA&I	<input checked="" type="checkbox"/>
DTIC TAB	<input type="checkbox"/>
Unannounced	<input type="checkbox"/>
Justification	
By	
Distribution/	
Availability Codes	
Dist	Avail and/or Special
A-1	



## LIST OF FIGURES

<u>Figure</u>	<u>Page</u>
1. Drawing of 11 Fibers to be Measured	2
2. Histogram of Area of 11 Fibers	15
3. Histogram of Perimeter of 11 Fibers	16
4. Histogram of Feret X of 11 Fibers	17
5. Histogram of DCircle of 11 Fibers	18

## LIST OF TABLES

<u>Table</u>	<u>Page</u>
1. Image Analysis Program for Measuring Fibers	13
2. Statistical Data (in Pixels) for 11 Fibers Measured	19

## PREFACE

In this report we describe the individual components and sequence of instructions of a computer program designed for the stereological measurement of fibers. Utilizing on-screen menu commands, we were able to design and execute a step-by-step analysis of a fiber measurement program. There was some difficulty in writing the program due to a lack of clear instructions in the operator's manual; however, conversations with company representatives helped us overcome these problems.



A COMPUTER PROGRAM FOR MEASURING FIBERS  
WITH THE ZEISS CSM 950 SCANNING ELECTRON MICROSCOPE

INTRODUCTION

The objective of this report is to describe a computer program designed for measuring stereological features of fibers by means of a Kontron image analysis system, which is built into a Zeiss CSM 950 Scanning Electron Microscope (SEM).

The program was written to permit interaction between the menu program displayed on the TV monitor and the keyboard or mouse. This increased capability permits the operator to change parameters within the measurement program depending on the type of fiber sample.

Since research presently being conducted in these laboratories on ballistic fibers, including biotechnologically derived fibers, requires information about fiber size, shape, etc., developing a menu-driven computer program facilitates that purpose.

METHODS

The sample model used was a drawing of 11 fibers (Fig. 1) which were visualized with a Dage MTI Model CCD 725 video camera with a Fuji C6 X 17.5 TV zoom lens. The images from the TV camera were transferred to a Zeiss CSM 950 SEM with a built-in Kontron image analysis system. The CSM 950 utilizes a CP/M operating system having a Z80 processor with a



FIGURE 1. Drawing of 11 Fibers to be Measured. The parameters measured (in pixels) were Area, Perimeter, Feret X (maximum projected length in the horizontal direction), and DCircle (area of fibers within a circle measuring 100 pixels in diameter).

20 Mbyte hard disk, 64 Kbyte random access memory (RAM) and a 16 Kbyte video RAM.

The image processing system performs evaluations in a series of steps. First, the image from a TV camera, SEM, light microscope or any other internal or external source is digitized and stored. Next, the grey image is processed for improved contrast and signal/noise ratio. Then, certain features are extracted, i.e., unwanted or overlapping features are eliminated, and features such as size, shape, etc., are analyzed. Finally, the derived data are processed statistically and a printout or photograph is obtained.

There are several major functions, which are grouped depending upon the specific sequence of a standard measuring procedure. The functions are:

1. Input - defines the type of input device and controls image digitization and storage.
2. Calibrate - controls the scaling factors, i.e., inches, microns, etc.
3. Enhance - rescales the grey levels, thereby enhancing image quality. Image noise reduction filters are used.
4. Image Edit - contains functions for the editing of images.
5. Segment - extracts the image background (distinguishes between objects and background).
6. Multiphase - enables further processing of binary (two grey levels) and multiphase (more than two grey levels) images.

7. Parameters - permits the selection of object-specific (single objects in a field of view), or field-specific (area of a visual field minus the area occupied by objects).

8. Evaluate - evaluates objects and outputs the results.

9. Advanced - performs geometric and arithmetic transformation of images, i.e., improves the quality of an image using Fast Fourier Transform methods.

10. Utilities - has general auxiliary functions.

11. Peripheral - permits control of peripheral devices, remote control, etc.

All the above functions have subfunctions that are displayed on the video monitor in the form of a menu. Dialogue with the program, allowing the operator to assign or delete numerical values as well as to control all of the system functions, is accomplished by utilizing the keyboard or a mouse.

The following is the selection sequence for the operator in order to obtain the desired program.

#### A. Calibrate

1. Scale - As all measurement data are indicated in the selected unit of calibration, the most important prerequisite for every geometric measurement is calibrating the measuring system. The calibration procedure is carried out according to the instructions in the operator's manual (Anonymous, 1986).

a. SCNO, the number of the scale factor to be activated, is set at 1.

## B. Evaluate

1. Reset - Before starting a new measuring sequence, e.g., after loading a new program, some conditions such as data buffering, identification number and scanning stage can be reset. The system then has the same status as after initialization.

a. Single (on) - Clears the buffer into which all the data from several fields are accumulated.

b. Class (on) - Clears the data summed up over several evaluated fields to provide average values (e.g., mean area percentage).

## C. Select Parameter

1. Area - Object specific area.

a. NCLS (20) - Number of classes when classification is selected. Classification refers to a range of measurement parameter values, i.e. area, perimeter, etc. into which classes (groups) of fibers would fit.

b. MODX (1-3) - Presentation and scaling mode of the X-axis in the histogram of the results.

c. MODY (1-10) - Presentation and scaling mode of the Y-axis in the histogram of the results.

d. Single (on) - Selection of a single list.

e. Class (on) - Selection of a classification.

f. Low 0.000 - Lower bound of classification.

g. High 5000 - Upper bound of classification.

## D. Select Parameter

1. Feret - Feret diameters (maximum projected lengths) in X and Y direction.

a. NCLS (20) - Number of classes when classification is selected.

- b. MODX (1-3) - Presentation and scaling mode of the X-axis in the histogram of the results.
- c. MODY (1-10) - Presentation and scaling mode of the Y-axis in the histogram of the results.
- d. Single (on) - Selection of a single list.
- e. Class (on) - Selection of a classification.
- f. Low (0.000) - Lower bound of classification.
- g. High (200.00) - Higher bound of classification.
- h. NCLY (20) - Number of classes when classification is selected.
- i. MDYX (1-3) - Presentation and scaling mode of the X-axis in the histogram of the results.
- j. MDYY (1-10) - Presentation and scaling mode of the Y-axis in the histogram of the results.
- k. LO.Y (0.000) - Lower bound of the classification.
- l. HI.Y (1000.0) - Upper bound of the classification.

#### E. Select Parameter

- 1. Perim - Object specific perimeter and length.
  - a. NCLS (20) - Number of classes when classification is selected.
  - b. MODX (1-3) - Presentation and scaling mode of the X-axis in the histogram of the results.
  - c. MODY (1-10) - Presentation and scaling mode of the Y-axis in the histogram of the results.
  - d. Single (on) - Selection of a single list.
  - e. Class (on) - Selection of a classification.

- f. Low (0.000) - Lower bound of classification.
- g. High (600.0) - Upper bound of classification.
- h. NCLT (10) - Number of classes when classification is selected.
- i. MDLX (1-3) - Presentation and scaling mode of the X-axis in the histogram of the results.
- j. MDLY (1-10) - Presentation and scaling mode of the Y-axis in the histogram of the results.
- k. IOLT (600.0) - Lower bound of classification.
- l. HILT (0.000) - Upper bound of classification.

#### F. Select Parameter

- 1. DCIRCL - Diameter of area-equivalent circle (area of fibers within a circle which measures 100 pixels in diameter).
  - a. NCLS (20) - Number of classes when classification is selected.
  - b. MODX (1-3) - Presentation and scaling mode of the X-axis in the histogram of the results.
  - c. MODY (1-10) - Presentation and scaling mode of the Y-axis in the histogram of the results.
  - d. Single (on) - Selection of a single list.
  - e. Class (on) - Selection of a classification.
  - f. Low (0.000) - Lower bound of classification.
  - g. High (100.00) - Upper bound of classification.

#### G. Select Utilities

- 1. LAB: Label definition - Labels are used as markers in the sequence of a measuring program to which a jump (analogous to a GOTO function) can be executed.
  - a. LBL#(1) - Identifier of the label, which serves as designation for a jump.

#### H. Input

1. TVON - Switches the color monitor directly to a TV signal.
  - a. Online (on) - Switches the online status on or off.

#### I. Select Utility

1. SYNC - The Image Processing System can be synchronized externally (e.g., by using a TV camera) or internally. In the internal mode, no TV input is possible.
  - a. Internal (off), Internal (on), or External (off).

#### J. Select Utility

1. Pause - Interrupts run of program.

#### K. Input

1. TVINP - This function stores:
  - a. TV Input - A TV image in the memory indicated by INP
  - b. INP (1) - Image memory into which the image is to be stored.

#### L. Choose Segment

1. DISC2L - This function separates objects from the background, by setting two thresholds. Either the grey levels inside or those outside the entered limits are set to grey value 0 (black) and constitute the background. Depending on the variable BINARY, the remaining object points either keep their original grey levels or are set to white (255).
  - a. INP (1) - Image memory containing the image to be processed.
  - b. OUT (2) - Image memory into which the discriminated image is stored.
  - c. LEV1 (0) - Lower discrimination limit.
  - d. LEV2 (118) - Upper discrimination limit.



- e. BINARY (on) - The object pixels become 255 (white), i.e. a binary image is produced.

#### M. Calibrate

1. MFRAME - This function is used to create a rectangular or circular measuring frame, which is sometimes necessary to correct errors (eliminate overlapping images) caused by the finite size of the field of view during the measurement process.

- a. INP (2) - Memory containing the image to be processed.
- b. WDSX (508) - X-side of the rectangular frame.
- c. WDSY (467) - Y-side of the rectangular frame.
- d. XO (4) - X position of the (rectangular or circular) frame.
- e. YO (6) - Y position of the (rectangular or circular) frame.
- f. RAD (100) - Radius of the circular frame.

#### N. Evaluate

1. IDENT - This function serves two purposes: Identification of discriminated images and elimination of objects following the conditions set by the measuring frame (function MFRAME, group Calibrate).

- a. INP (2) - Image memory containing the image to be identified.
- b. OUT (3) - Memory containing the identified and frame-corrected image.
- c. MARG (1) - Controls the meaning of the measuring frame defined by MFRAME (function group Calibrate). This determines whether or not the fiber touching the rectangular or circular frame on the video monitor is measured. For example, one might not want to measure a fiber extending beyond the frame's margin.

d. 8-CONN (on) - Method of identification, on = 8 CONN.

Q. Evaluate

1. MEASURE - This function executes the measurement of selected parameters such as area, perimeter, etc.

- a. OBJ - Must be switched on if object-specific parameters are among the previously selected parameters (AREA-TIME).
- b. INP (3) - Image memory containing the discriminated OBJ image.
- c. GRIM (1) - Image memory containing the original grey image.
- d. AUX1 (9) - Auxiliary image memory.
- e. AUX2 (9) - Auxiliary image memory.
- f. CHAN (1) - Channel selection.
- g. SPEC (1) - Distance in pixels between consecutive lines of the "line grid". The lines are simply lines of the image memory, or pixel addresses with varying X, but constant Y coordinates.
- h. FLG# (1) - Logical flag indicating a buffer overflow.
- i. REJECT (on) - Interactive rejection of undesired objects, preceding an automatic measurement.

P. Evaluate

1. OUTCLS - (Output of Classification) This function displays the results of the classification in the form of a list and/or as a histogram.

- a. HISTO (on) - Provides measured data in histogram form.
- b. HALT (on) - Presents data on the monitor; the sum histogram can be manipulated.

c. NFAC (0.000) - Defines the area (in user units) into which the absolute counts should be normalized. A value of 0.0 gives normalized absolute frequencies.

Q. Evaluate

1. OUTSGL - (Output of single list) This function displays the individual lists of all measured data.

a. HALT (on) - Displays lists on the data monitor after each measurement.

R. Select Utility

1. Pause - This function interrupts the running of a measuring program.

## RESULTS

The computer program designed and written for the stereological analysis of several fibers (Table 1) was utilized to measure their Area, Perimeter, Feret X and D Circle. The inputs to the parameters of this program were as follows:

Area (high upper bound of classification) - 5000

Feret X (higher bound of classification) - 200

Feret X (upper bound of classification) - 1000

Perimeter (high upper bound of classification) - 600

DCircle (high upper bound of classification) - 100

The printout of the computer program can be seen in Table 1, the plots for each parameter can be seen in Figures 2 to 5 and the printout of the statistical data can be seen in Table 2.

## CONCLUSIONS

A computer program was designed and written to measure certain stereological parameters of fibers and to generate subsequent statistical analysis of the results. The program was developed to be driven by a series of menu commands, which allowed for direct interaction between the image appearing on the video monitor and the keyboard.

TABLE 1. Image Analysis Program for Measuring Fibers Using the Zeiss CSM 950 SEM with a Built-in Kontron Image Analyzer.

MAJOR FUNCTIONS	SUB FUNCTIONS	INTEGER VALUES	VARIABLES	SCALE VALUES
A. CALIBRATE	SCALE	SCNO 1		
B. EVALUATE	RESET			
			SINGLE CLASS	
			PARAM	
C. PARAMETERS	AREA	NCLS 20 MODX 1 MODY 1	SINGLE CLASS	LOW 0.000 HIGH 5000.0
D. PARAMETERS	FERET	NCLS 20 MODX 1 MODY 1 NCLY 20 MDYX 1 MDYY 1	SINGLE CLASS	LOW 0.000 HIGH 200.0 LO.Y 0.000 HI.Y 1000.0
E. PARAMETERS	PERIM	NCLS 20 MODX 1 MODY 1 NCLT 10 MDLX 1 MDLY 1	SINGLE CLASS	LOW 0.000 HIGH 600.0 LOLT 600.0 LILT 0.000
F. PARAMETERS	DCIRCL	NCLS 20 MODX 1 MODY 1	SINGLE CLASS	LOW 0.000 HIGH 100.0
G. UTILITIES	LAB1	LBL# 1		
H. INPUT	TVON	ONLINE		
I. UTILITIES	SYNC			

(Continued)

TABLE 1. Image Analysis Program for Measuring Fibers Using the Zeiss CSM 950 SEM with a Built-in Kontron Image Analyzer. (Continued)

MAJOR FUNCTIONS	SUB FUNCTIONS	INTEGER VALUES	VARIABLES	SCALE VALUES
J. UTILITIES	PAUSE			
K. INPUT	TVINP	INP 1		
L. SEGMENT	DISC2L	LEV1 0 LEV2 118	BINARY LEV1 0	
M. CALIBRATE	MFRAME	INP 2 WDSX 508 WDSY 467 XO 4 YO 6 RAD 100		
N. EVALUATE	IDENT	INP 2 OUT 3 MARG 1	8-CONN	
O. EVALUATE	MEASUR	INP 3 GRIM 1 AUX1 9 AUX2 9 CHAN 1 SPAC 1 FLG# 1	OBJ  REJECT	
P. EVALUATE	OUTCLS		HIST	NFAC 0.000
			HALT	
Q. EVALUATE	OUTSGL		HALT	
R. UTILITIES	PAUSE			
TOTAL PROGRAM LENGTH		340 BYTES		

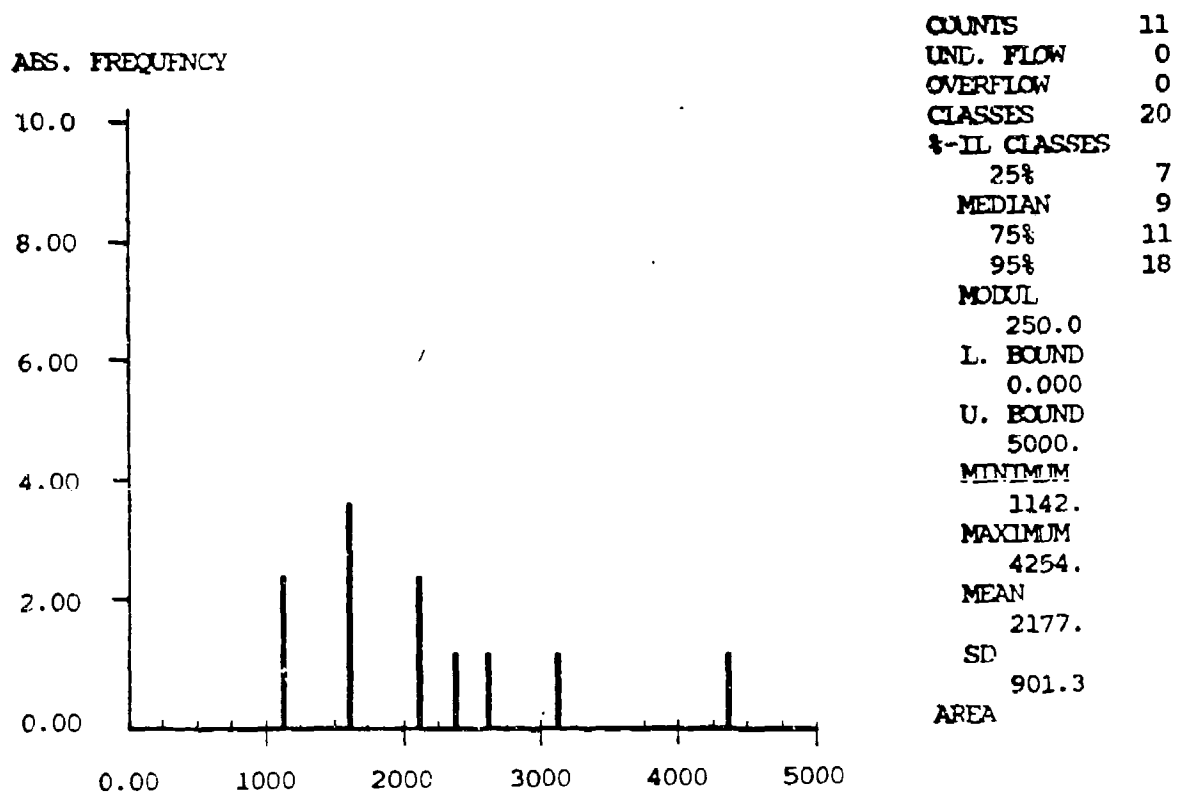


FIGURE 2. Histogram of Area Measurements from Figure 1.  
 The X axis values are in pixels.  
 The Y axis represents the frequency distribution.

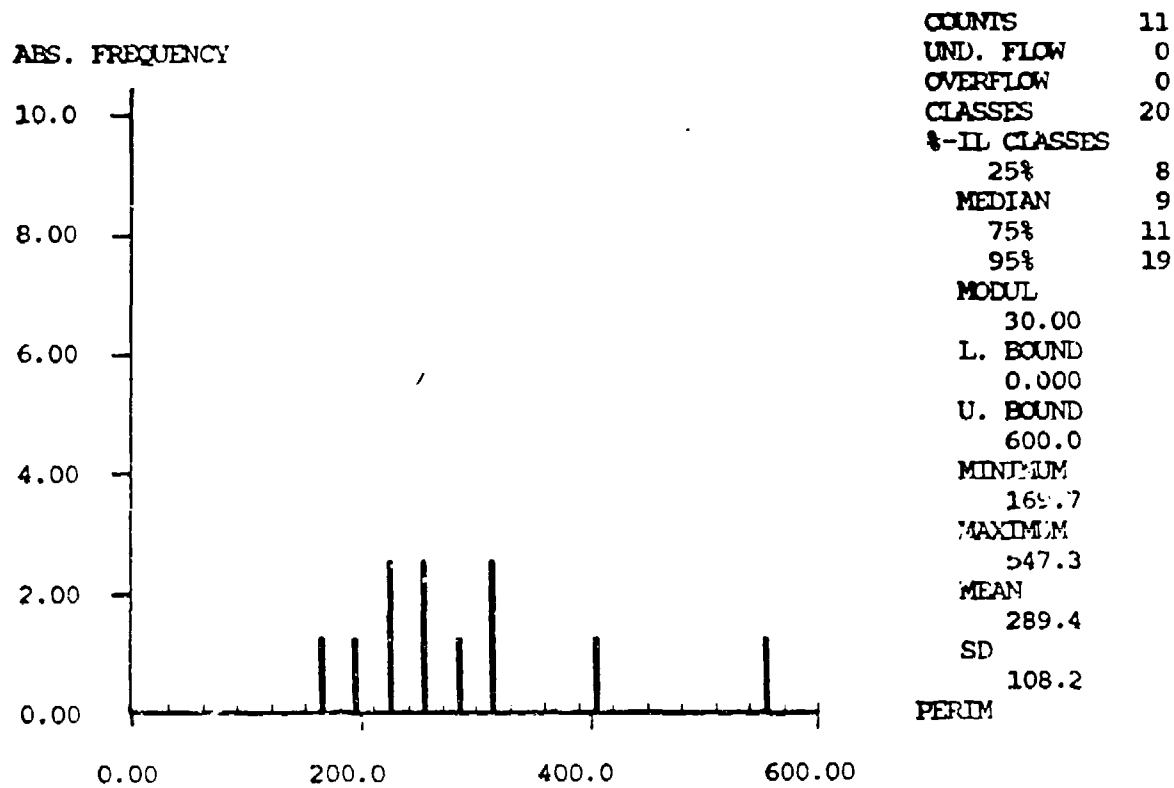


FIGURE 3. Histogram of Perimeter Measurements from Figure 1.  
 The X axis values are in pixels.  
 The Y axis represents the frequency distribution.



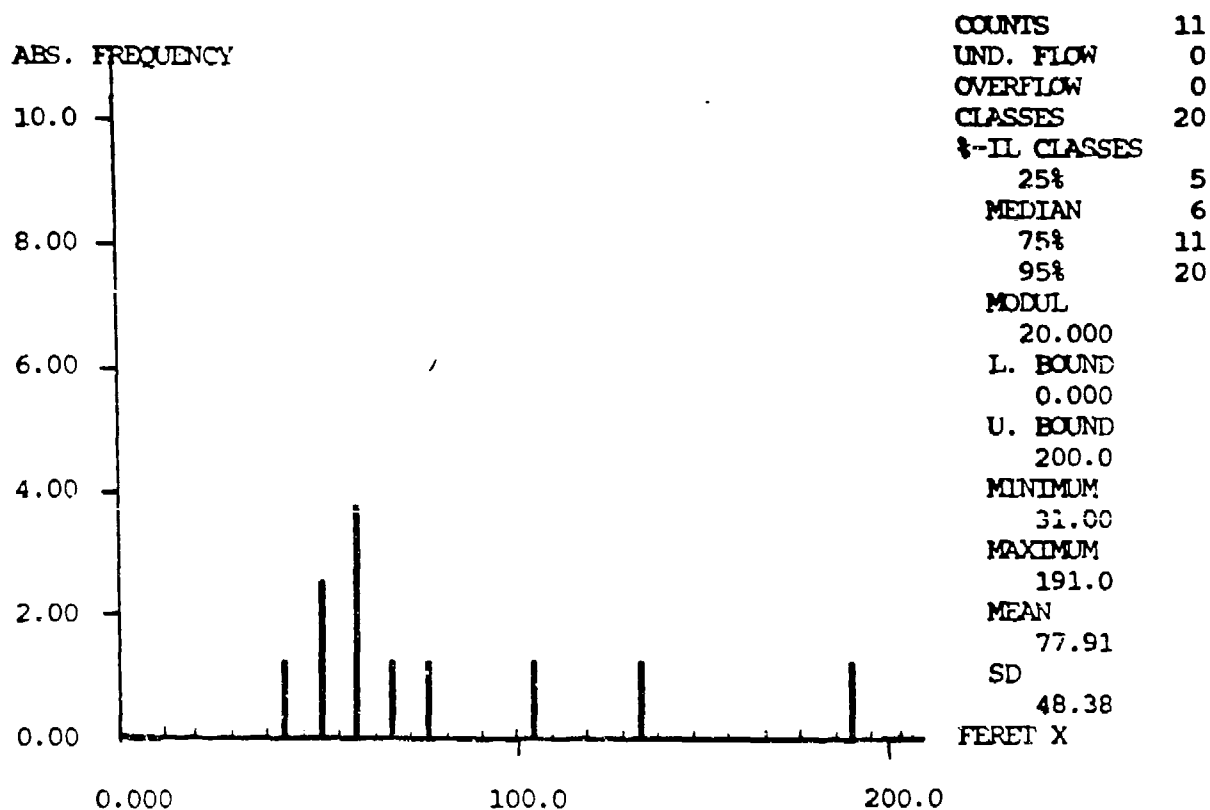


FIGURE 4. Histogram of Feret X Measurements from Figure 1.  
 The X axis values are in pixels.  
 The Y axis represents the frequency distribution.

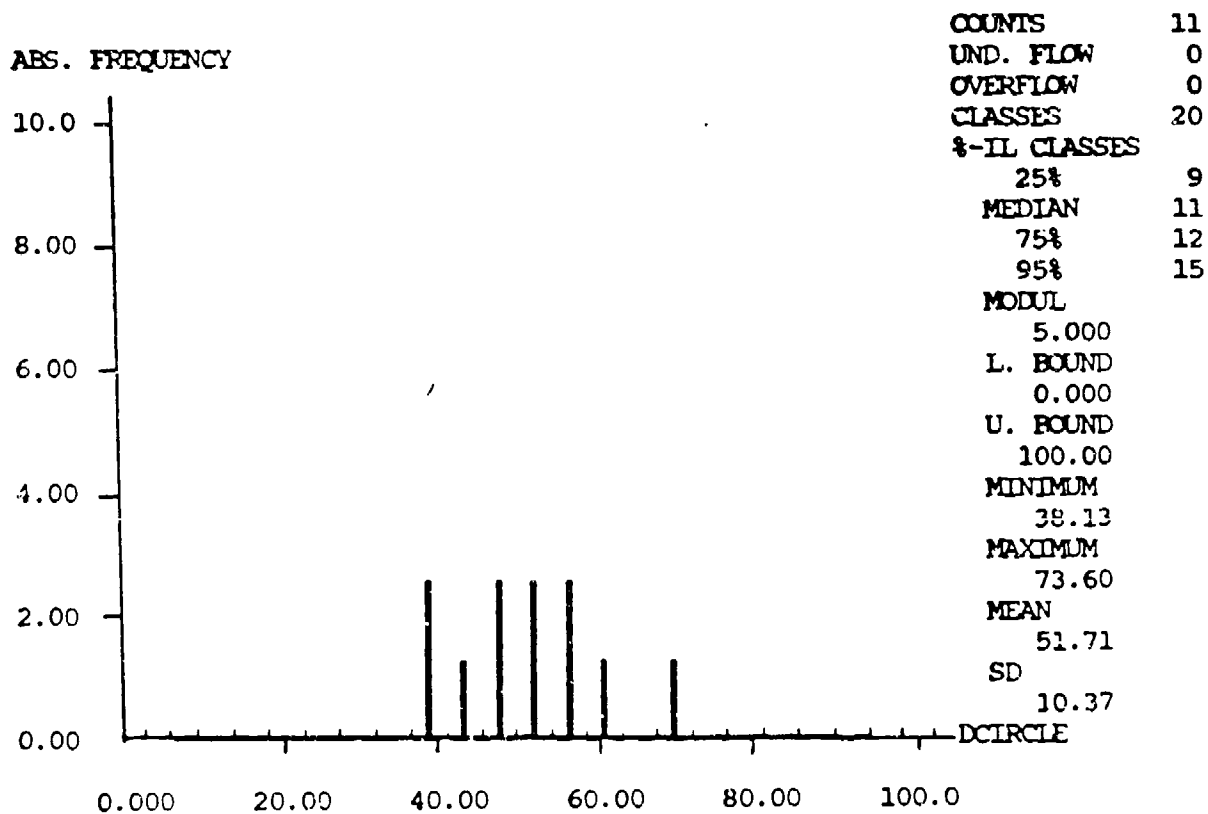


FIGURE 5. Histogram of an Area Equivalent to a Circle (called a DCircle) within Figure 1.  
The X axis values are in pixels.  
The Y axis represents the frequency distribution.

TABLE 2. Statistical Data (in Pixels) for 11 Fibers Measured

<u>COUNT</u>	<u>AREA</u>	<u>PERIM.</u>	<u>FERETX</u>	<u>DCIRCLE</u>
1	3111.0	413.4	137.0	62.94
2	2134.0	270.0	53.0	52.13
3	2512.0	318.4	58.0	56.55
4	1731.0	250.8	103.0	46.95
5	2378.0	302.6	52.0	55.03
6	1691.0	226.0	41.0	46.40
7	2196.0	269.5	44.0	52.88
8	1142.0	169.7	68.0	38.13
9	1217.0	196.7	79.0	39.36
10	1582.0	219.2	31.0	44.88
11	4254.0	547.3	191.0	73.60

# REFERENCE

Anonymous, 1986. "Kontron SEM-IPS Operator's manual," Volume II, Release 4.4, 410 pages.

This document reports research undertaken at the US Army Natick Research, Development and Engineering Center and has been assigned No. NATICK/TR-91/DC4 in the series of reports approved for publication.